WORK HOOK INSTRUCTIONS FOR USE

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1. INTRODUCTION

This document contains indications for use and assembly of the work hook. Please read it carefully before using the product to ensure that you are aware of the safety instructions and to inform your patients.

The ORTHOPUS Work Hook is to be used solely for the exoprosthetic fitting of the upper extremities in combination with other arm components.

2. INDICATIONS FOR USE

The ORTHOPUS Work Hook is designed to be part of an upper-limb body powered prosthetic.

Rock solid, functional, easy to use and maintain, this cable-operated Work Hook makes a great variety of manual work possible.

The ORTHOPUS Work Hook is available in the following options:

- Left Work Hook with a M12x1.75 thread (reference: OR-011201L)
- Right Work Hook with a M12x1.75 thread (reference: OR-011201R)
- Left Work Hook with a ½"-20 UNF thread (reference: OR-011202L)
- Right Work Hook with a ½"-20 UNF thread (reference: OR-011202R)
- Left Work Hook with a M10x1.5 thread (reference: OR-011203L)
- Right Work Hook with a M10x1.5 thread (reference: OR-011203R)

The grip force of the hook can be adjusted by changing the number of elastic bands. The more elastic bands are added, greater will be the grip force. However, if the grip force is high the user will have to use a lot of energy to open the hook and it could lead to fatigue.

Few examples of object grips are shown on Page 2 to illustrate some use cases of the hook.

3. SAFETY INSTRUCTIONS

Please make your patients aware of the following safety instructions:

Warnings :

- Do not expose to excessive temperatures or forces
- Do not use for extreme or unusual activities (such as extreme sports)

- Users must comply with local regulations on the operation for motorized vehicles in his/her country
- Do not expose the prosthesis or his components to liquids, dust, acids, aggressive cleaning agents or flames

Precautions for use :

- In case of unusual stress or shock, the prosthesis must be inspected by a prosthetist
- Risk of injuries due to improper use

4. TECHNICAL SPECIFICATIONS

Work Hook	
Weight (g)	112 g
Total height (mm)	140 mm
Max opening width (mm)	120 mm
Pinch force (40 mm opening width, 1 elastic band)	1.7 kg (16.7 N)
Grip force (80 mm opening width, 1 elastic band)	1.6 kg (15.7 N)

5. ASSEMBLY

As the Work Hook is activated by a cable, a cable with proper termination is required. This one needs to be inserted at the end of the Work Hook Thumb (dimensions given on Fig. 2).

The Work Hook is compatible with the ORTHOPUS Universal Wrist by screwing it into the tool holder (reference: OR-0109AM0x). In any case you wish to mount the Work Hook on a tool holder, add some glue on the Work Hook thread as shown on Fig. 3.

With the three different threads proposed, M12x1.75 (reference: OR-011201x), ½"-20 UNF (OR-011202x) or M10x1.5 (OR-011203x), the Work Hook is compatible with most of the standard components of the market.

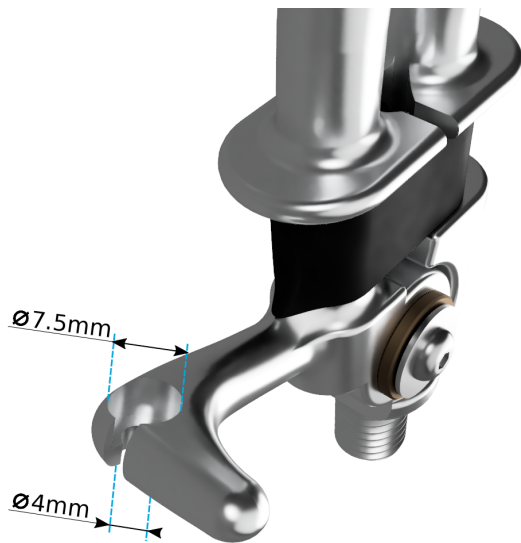
Figure 2

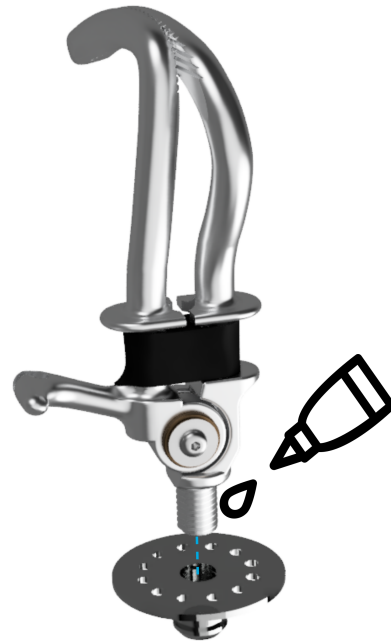
Figure 3

6. HOOK GRIP

In order to improve grasping capacities of the Work Hook, it is possible to add a silicone grip at both extremities of the hook.

To add the silicone grips, insert the Work Hook extremities inside the grips. If the insertion is too difficult, you can apply some soapy water on the hook. At the end of the insertion process, the equipped Work Hook should look like the one presented here:

With the silicone grips, it might be difficult to grasp small and fine objects (a sewing needle for example) with the tip of the Work Hook. If the fine grasp function is essential to the user, silicone grips can be cut in order the hook can be fully closed (as shown below). Another solution is to not use the silicone grips.

7. MAINTENANCE

If cleaning or unjamming of the Work Hook is required, it can be fully disassembled by unscrewing the screw 3 (see Fig. 4). As glue was used to lock the screw in the hook, it might be difficult to unscrew it. Please be careful to use appropriate tools and not to damage the screw head. To ensure a robust screw mounting, add a little drop of glue in the screw hole before tightening it.

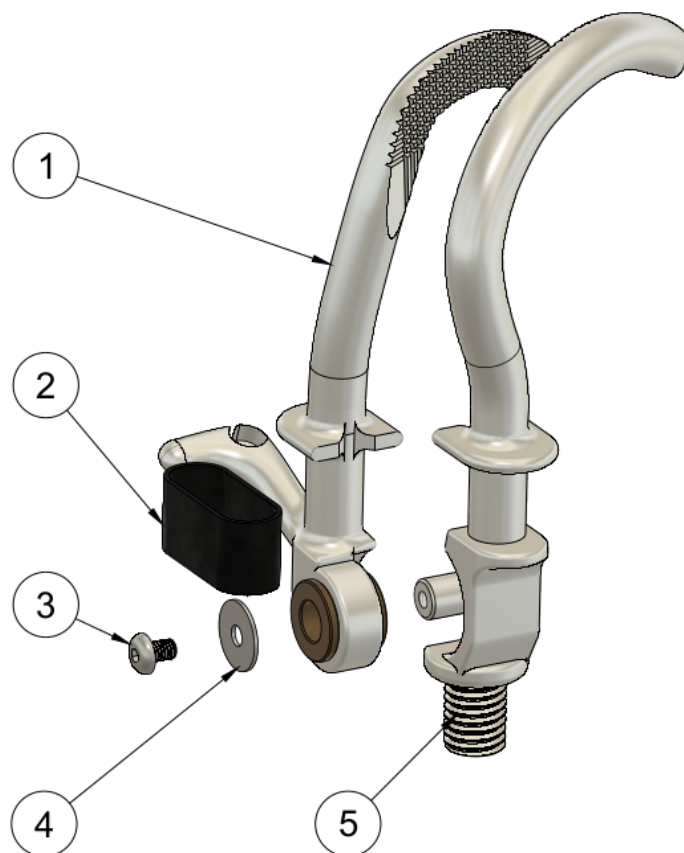
Figure 4

8. REGULATORY INFORMATION

All legal conditions are subject to the respective national laws of the country of use and may vary accordingly.

Responsability:

As manufacturer, ORTHOPUS will not assume liability for damage caused by disregard of this document, particularly due to improper use or unauthorised modification of the product.

CE Marking:

This product meets the requirements of Regulation (EU) 2017/745 for medical devices. This product has been classified as a class I device according to the classification criteria outlined in Annex VIII of the directive.

9. DISPOSAL

This product, its components and its packaging should be disposed of in accordance with respective national environmental regulations. Observe the information provided by the responsible authorities in your country regarding return, collection and disposal procedures.

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